


BRODA Tranquille Glider



A gentle gliding action creates a calming, soothing effect for agitated individuals. Particularly useful for Alzheimers and dementia patients as the gliding motion provides stimulation for active patients and helps reduce the likelihood of wandering.

The auto-locking glider incorporates a falls prevention mechanism in that it only begins to glide once the resident sits in the chair and automatically stops as the chair is vacated. It also features a caregiver lock to prevent unwanted gliding.

The Glider also features the Broda Comfort Tension Seating which provides long term seating comfort as well as pressure redistribution and accommodates conditions such as kyphosis and lordosis.

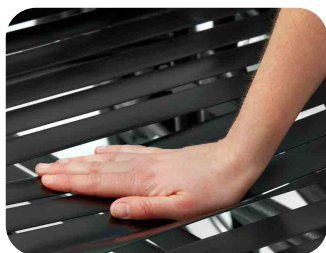

The Unique Comfort Tension Seating system provides therapeutic pressure relief which aids pressure ulcer treatment and pain management.

Broda Comfort Tension Seating

Powder coated tubular steel frame

Height adjustable armrests

Auto-locking mechanism

- falls prevention

Caregiver lock

- makes chair stationary

Lumbar support

Side panels, seat pad and back pad



Seat Width	Seat Depth	Seat Height	Back Height	Armrest Height	S.W.L.
51cm	43cm	46cm	66cm (Low-back) 81cm (High-back)	19 - 34cm	113kg