

Malibu Care Bed

Ideal for use in Aged Care facilities, the Malibu Care Bed is also a perfect option for using in home environments. The Malibu covers all the necessary functionality of a clinical bed, with standard features included that have taken the needs of users and carers into consideration. Hoist friendly, a Safe Working Load of 220kg and convenient pedal central-locking on all four twin wheel castors.




Available in **Single** and **King Single**.



Rotating wall protection **buffers** on head end of bed.

Pedal **central-locking** on all four twin wheel castors.

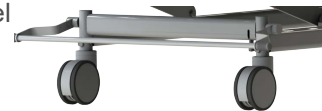

Easy to use **handset control** adjustment using **Linak electronic** drive system.

$\frac{3}{4}$ Tubular fold-down horizontal **handrails** can be easily removed or clamped on.



Mattress retainers on four corners prevent your mattress from shifting when bed is adjusted.



The built-in bed **extension** adds an extra 19.5cm to the mattress base length.



The Malibu comes with a 4-Sectioned mattress platform. This allows for auto contouring as well as:

Independent adjustment of

- Mattress platform height (28.5 - 68.5cm)
- Backrest angle (0° - 70°)
- Knee break angle (0° - 35°)


The sheet steel deck panels are easily removable to facilitate maintenance and servicing.



The head and footboards on the Malibu Care Bed come standard with a stylish neutral beech timber finish.



OPTIONAL ACCESSORIES

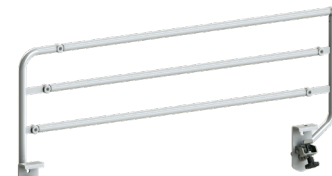

¾ Tubular horizontal rails



Battery back-up unit



Self-help pole



Self-Assist rail - fold down



IV pole

	Width - overall Width - platform	Length - overall Length - platform	Extension	Height range	Backrest angle	Knee break angle	S.W.L.
Single	92cm 88cm	214cm 202cm	19.5cm	28.5 - 68.5cm	0° - 70°	0° - 35°	220kg
King Single	109cm 107cm	214cm 202cm	19.5cm	28.5 - 68.5cm	0° - 70°	0° - 35°	220kg